

Refraction unit X2S



X2S



Refraction unit X2S

- Sliding table for two instruments
- Adjustable height of table and chair
- Column for chart projector with adjustable LED overhead lamp
- Built-in control panel with ambient light
- Wheelchair friendly
- Drawer for trial lenses set, swiveling and height adjustable
- Easy right-/left- handed transformable at customer's place
- Adjustable adapter for chin rest

Customization options

- Color choices: The X2S unit is available in a wide range of colors for its metal parts, tabletops, and chair upholstery
- Tabletop customization: The size and shape of the tabletops can also be tailored to suit the specific needs and preferences

Optional

- Foot rests & arm rests: For added comfort during extended examinations
- Phoropter arms: can be equipped for a comprehensive setup
- Monitor holder: For convenient display and positioning of monitors during patient consultations.
- Foot switch control: For adjusting the table and chair heights conveniently during examinations.
- Factory-Installed data cables: For seamless instrument integration with FireWire, USB, or RS232 connectivity.
- Reclining patient chair: manual or motorized
- Patient chair sliding assembly: manual or motorized

X2S

1. Chair with reclining back.



3. Built-in control panel with ambient light.





2. Adjustable adapter for chin rest.



Technical information

Table height from the ground	between 70 cm to 100 cm +/-2mm
Max tabletop load capacity	60 kg (1st position – max 20kg, 2nd position – max 40kg)
Max chair load capacity	150 kg
Instrument power supply	ranges from 3.5V to 24V AC&DC, adaptable to different instrument power requirements
Power supply	230V, 50/60Hz, max 1.4A. 115V version available on request
Weight*	aprox 150kg
Max weight of the instruments (1st – 2nd)	15kg – 30kg

*Varies depending on options

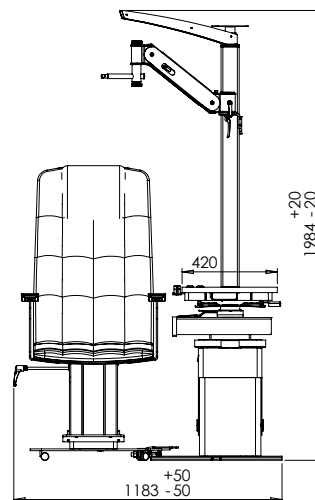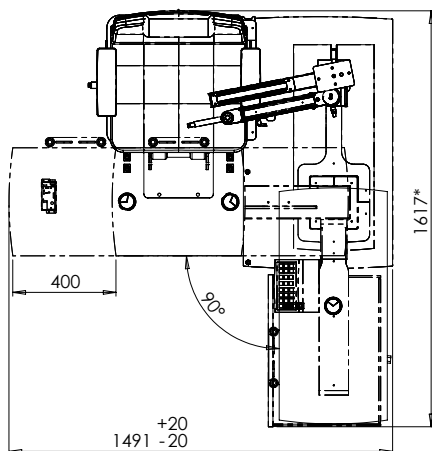

*Varies depending on options



WWW.MEDINSTRUS.COM
INFO@MEDINSTRUS.LT

medinstrus